



## Table of Contents

Stainless Steel Tables	05
Mayo Stands	13
Additional OR Equipment	19
Storage, Transfer and Transport	23
Accessories for OR Tables	29
Arm and Hand Surgery Tables	39
IV Poles and Accessories	47
Index	56

# Stainless Steel Tables

A flat surface impervious to infection supported by 2 or 4 legs, used to hold instruments in the operating room.



# Stainless Steel Tables



**sayno2rust.**<sup>®</sup>

Look for this logo to  
see eligible products



sayno2rust Caster



sayno2rust Caster  
with Brake

## REVEALED EDGE

Our tables have an 11/2" revealed edge. This greatly enhances the rigidity of the top and shelf. The more material you turn down, the stronger the edge will be. What does that mean for you? *The edge of your table is stronger so it will stand up to use in the OR better.*

## MATERIAL THICKNESS

Mid Central Medical Tables are 16 gauge construction. That means the top and shelf are constructed of 16 gauge (0.0625") material instead of 18 gauge (0.05"). That is 25% thicker material! What does that mean for you? *Less deflection of the top when you place heavy items on it.*

## REINFORCING SUPPORTS

Welded in reinforcing supports under the top and shelf of any table longer than 4' with sound dampening cushion to prevent vibration of the top during transport.

## SQUARE LEGS

Many customers ask, why square instead of round? *The square leg has more surface for welding the shelf and top to the leg. 1" Square tubing also has 25% more actual stainless steel material in it than 1" round tube.*

## CASTER UPGRADE

Are you tired of seeing rust in your Operating Room? We can help! Now all our instrument tables can be upgraded with our trademarked **sayno2rust** casters. These are durable, easy to clean casters with soft-gripped non-marking tread. Available with brake and all swivel wheels. Simply add the item numbers below as a line item on your purchase order and **sayno2rust!**

CSTR-TBLE	sayno2rust Composite Caster Upgrade, 4 swivel
CSTR-TBLEBRK	sayno2rust Composite Caster with Brake Upgrade, 2 braking, 2 swivel

# Stainless Steel Tables



## HCM500 Series

### Instrument Tables with Shelf

Our most common instrument table with lower shelf. Also known as a back table in the OR. Fully welded stainless steel construction. These tables are manufactured to our OR specifications and do not have any open cavities to collect unwanted substances. All edges are glove friendly. Ships fully assembled.

- › New expanded size range
- › No ripped gloves
- › Caster upgrade available

**sayno2rust.**

HCM501	16"D x 20"W x 34"H
HCM502	16"D x 30"W x 34"H
HCM503	18"D x 33"W x 34"H
HCM519	20"D x 20"W x 34"H
HCM504	20"D x 36"W x 34"H
HCM506	20"D x 48"W x 34"H
HCM551	24"D x 30"W x 34"H
HCM505	24"D x 36"W x 34"H
HCM507	24"D x 48"W x 34"H
HCM508	24"D x 60"W x 34"H
HCM509	24"D x 72"W x 34"H
HCM303634SHCA	30"D x 36"W x 34"H
HCM304834SHCA	30"D x 48"W x 34"H
HCM306034SHCA	30"D x 60"W x 34"H
HCM307234SHCA	30"D x 72"W x 34"H

\*Custom sizes available



## HCM510 Series

### Instrument Tables with H-brace

Instrument table with H-brace adds more clearance around the base of the tables for storage and fitting around obstacles. Fully welded stainless steel construction. These tables are manufactured to our OR specifications and do not have any open cavities to collect unwanted substances. All edges are glove friendly. Ships fully assembled.

- › New expanded size range
- › No ripped gloves
- › Caster upgrade available

**sayno2rust.**

HCM510	16"D x 20"W x 34"H
HCM511	16"D x 30"W x 34"H
HCM512	18"D x 33"W x 34"H
HCM519-H	20"D x 20"W x 34"H
HCM513	20"D x 36"W x 34"H
HCM515	20"D x 48"W x 34"H
HCM551-H	24"D x 30"W x 34"H
HCM514	24"D x 36"W x 34"H
HCM516	24"D x 48"W x 34"H
HCM517	24"D x 60"W x 34"H
HCM518	24"D x 72"W x 34"H
HCM303634HBCA	30"D x 36"W x 34"H
HCM304834HBCA	30"D x 48"W x 34"H
HCM306034HBCA	30"D x 60"W x 34"H
HCM307234HBCA	30"D x 72"W x 34"H

\*Custom sizes available



## HCM560 Series

### Nested Instrument Tables

Instrument tables that nest together for storage. U-brace also allows for maneuverability around obstacles. Fully welded stainless steel construction. These tables are manufactured to our OR specifications and do not have any open cavities to collect unwanted substances. All edges are glove friendly. Ships fully assembled.

- › Less storage space required
- › No ripped gloves
- › Caster upgrade available

**sayno2rust.**

HCM560	Set of 6 Nested Tables
HCM561	14"W x 28)L x 32"H
HCM562	16"W x 32)L x 34"H
HCM563	18"W x 36)L x 36"H
HCM564	20"W x 40)L x 38"H
HCM565	22"W x 44)L x 40"H
HCM566	24"W x 48)L x 42"H

\*Custom sizes available

# Stainless Steel Tables



## HCM530 Series

### Height Adjustable Instrument Table

Tables designed to straddle the OR Table for use during cases. They adjust in height from 36" to 56" to accommodate the height of the OR Table. They are fully welded construction, with no exposed bolts or protrusions on the top to puncture the drape. Handle removes and stores inside the lower leg if needed. Can be operated from either end. All exposed surfaces are of stainless steel construction. Reinforcing channel under top to provide extremely good rigidity. Top material is 14 gauge. Rolls on 4" premium casters, 2 have top brake to secure in place. Can also just be used as an ergonomic back table that can be adjusted in height.

- › **No more punctured drapes**
- › **Easy to clean**
- › **Ease of mobility and locking in place**
- › **Ergonomic for reduced body fatigue**

HCM530	20"W x 36"L x 36-56"H
HCM531	24"W x 36"L x 36-56"H
HCM532	24" W x 48" L x 36-56" H
HCM533	46" W x 34" L x 36-56" H
HCM535	24" W x 60" L x 36-56" H

\*Custom Sizes available



## HCM550 Series

### Two-tier Back Table

This two tier back table makes a great addition to any surgical suite, especially where larger cases are performed and floor space is tight. Constructed of 16 gauge stainless steel, with 1 1/2" square legs makes this an extremely sturdy table. Top shelf has 16" clearance to top of table. The 10 degree tilt on the upper shelf makes it easy to view instruments that are on the upper level.

- › **Available in 3 different sizes**
- › **4" locking casters on front 2 legs**
- › **30" depth allows for extra placement of instruments**

HCM554	30"D x 48"W x 34"H
HCM555	30"D x 60"W x 34"H
HCM556	30"D x 72"W x 34"H

\*Custom Sizes available



## HCM500 Series

### Drop-Leaf Instrument Tables

Stainless Steel Instrument table with a 12" extension on either side. With these drop down extensions, you essentially get 3 different tables in the same unit! Whether it is a small case or large, you can adjust the table size to suit your needs.

- › **Less space to store**
- › **Fit to the size case you are performing**
- › **No expensive custom drapes needed**
- › **Caster upgrade available**

**sayno2rust®**

HCM505-DPL	24"W x 36-60"L x 34"H
HCM507-DPL	24"W x 48-72" L x 34" H

# Stainless Steel Tables

## HCM520 Series

### Utility Tables

Utility Table with 3-sided guardrail and lower drawers. Standard size is 16"W x 20"L and 34"H, custom sizes available. 2 locking casters on front to secure unit in place. Locking top drawer available. Fully welded construction and made to OR specifications.

- › Drawers slide freely
- › No ripped gloves
- › Custom options available, get what you want
- › Caster upgrade available

**sayno2rust®**

HCM520	One Drawer and Lower Shelf
HCM521	Two Drawers and Lower Shelf
HCM522	Three Drawer and Lower Shelf
HCM523	Four Drawers and Lower Shelf



If you want to make sure something doesn't fall off, we can add a **guardrail** to any table.



Need more storage? A drawer or **drawers** can be added to any table.



3" **sayno2rust** composite caster upgrade.

# Stainless Steel Tables



## HCM523-T Series Stainless Steel Anesthesia Cart

HCM523-T

Stainless Steel  
Anesthesia Cart

This Anesthesia Cart boasts the largest tray size to hold instruments during spinal cases. Four lower drawers as well as a lower shelf provide extra storage. 2 locking casters to secure when needed.

- › Fully welded stainless steel construction
- › Large removable tray
- › Lots of storage space
- › Caster upgrade available

**sayno2rust®**



## HCM550 Series Stainless Steel Utility/Prep Table

HCM550

16"W x 20"L x 34"H

\*Custom sizes available

Utility Table with 4-sided guardrail around top and lower shelf. Standard size is 16"W x 20"L and 34"H. Also available in a larger size; custom sizes and configurations available on request. 2 locking casters on front to secure unit in place. Fully welded construction and manufactured to our medical specifications; suitable for use in the OR.

- › Heavy Duty 16 gauge #304 stainless steel construction
- › 4-sided guardrail
- › No open cavities to collect unwanted substances
- › Caster upgrade available

**sayno2rust®**



## HCM520 Series 18" x 33" Utility/Prep Table

HCM526

One Drawer and Shelf

HCM527

Two Drawers and Shelf

HCM528

Three Drawers and Shelf

Utility Table with 3-sided guardrail and lower drawers. Standard size is 18"W x 33"L and 34"H. Lower shelf adds extra storage space. 2 locking casters on front to secure unit in place. Locking drawers available. Fully welded construction and manufactured to our medical specifications; suitable for use in the OR.

- › Heavy Duty 16 gauge #304 stainless steel construction.
- › Drawers mounted on ball bearing glides.
- › 3-sided guardrail.
- › No open cavities to collect unwanted substances
- › Caster upgrade available

**sayno2rust®**

# Stainless Steel Tables



## HCM520 Series

### 20" x 38" Utility/Prep Tables

HCM524 Two Drawers and Shelf

HCM525 Four Drawers and Shelf

Utility Table with 3-sided guardrail and lower drawers. Standard size is 20"W x 38"L and 34"H. Also available in a larger size; custom sizes and configurations available on request. 2 locking casters on front to secure unit in place. Locking drawers available. Fully welded construction and manufactured to our medical specifications; suitable for use in the OR.

- › **Heavy Duty 16 gauge #304 stainless steel construction.**
- › **Drawers mounted on ball bearing glides.**
- › **3-sided guardrail.**
- › **No open cavities to collect unwanted substances**
- › **Caster upgrade available**

**sayno2rust.**



## HCM540 Series

### Instrument Tables with Drawers

HCM548 24"W x 48)L x 34"H, Two drawers

HCM549 24"W x 72)L x 34"H, Three drawers

Stainless Steel Instrument Tables with added storage. Either 2 or 3 drawers under top of table. Lower shelf adds extra storage space. Fully welded stainless steel construction.

- › **2 locking casters, 2 swivel**
- › **No ripped gloves**
- › **Caster upgrade available**

**sayno2rust.**



## HCM590 Series

### Central Sterile Instrument Wrapping Station

HCM596 32"D x 60)L x 34"H

This stainless steel instrument wrapping station has extra large plastic bins. With the stainless steel top and rack for wrapping instruments in sterile wrap, this table is everything you need. Ergonomic 35" work height and heavy duty 4" caster that lock in place ensure this table will last for many years. Standard size is 32" x 60" x 35"high. Can be easily adjusted in size to meet your needs.

- › **Sturdy work style table design.**
- › **Removable bins for easy cleaning and reloading instruments.**
- › **Easy to move but can still lock in place with locking casters**

# Stainless Steel Tables

## HCM500 Series

### Work Tables

These work tables are designed for working on with much more strength than standard instrument tables. A great addition to any OR, Central Sterile Supply, Laboratory or Pharmacy. The top is comprised of 14 gauge material and a reinforcing channel on the underside to prevent deflection. Supported by 1.5" 16 gauge stainless steel square legs make it very sturdy. Table can be configured with either leg levelers or 4" casters that allow the table to be locked in place during use. Ships fully assembled and ready to use.

- › **Made exclusively for the medical industry with glove friendly edges and no cavities to harbor bacteria**
- › **2" revealed edge with return to create a very sturdy table edge**
- › **Available with:**
  - 4" caster, two swivel locking, two swivel
  - Leg levelers that adjust in height to compensate for uneven floor substrates
  - With lower shelf for additional storage
  - With H-brace for extra clearance around base



Leg Leveler



H-Brace

HCM570	16"W x 20"L x 35"H
HCM571	16"W x 30"L x 35"H
HCM572	18"W x 33"L x 35"H
HCM573	20"W x 36"L x 35"H
HCM579	20"W x 48"L x 35"H
HCM574	24"W x 36"L x 35"H
HCM580	24"W x 48"L x 35"H
HCM582	24"W x 60"L x 35"H
HCM583	24"W x 72"L x 35"H
HCM589	24"W x 96"L x 35"H
HCM578	30"W x 30"L x 35"H
HCM588	30"W x 36"L x 35"H
HCM581	30"W x 48"L x 35"H
HCM584	30"W x 60"L x 35"H
HCM586	30"W x 72"L x 35"H
HCM590	30"W x 84"L x 35"H
HCM591	30"W x 96"L x 35"H
HCM592	36"W x 36"L x 35"H
HCM593	36"W x 48"L x 35"H
HCM585	36"W x 60"L x 35"H
HCM587	36"W x 72"L x 35"H
HCM594	36"W x 84"L x 35"H
HCM595	36"W x 96"L x 35"H

Customize your order: Item # + option

With Shelf add "-S"

With H-brace add "-H"

With Casters add "-CA"

With Leg Levelers add "-LL"



## Mayo Stands

Mobile Instrument Stand with removable tray that provides a place for sterile instruments and supplies used during surgery that is positioned over or adjacent to the surgical site.

# Mayo Stands

## HCM750 Series

### Foot Control Mayo Stands

Our most common Mayo Stand. Time tested and proven design. Removable tray for easy cleaning. Stainless steel foot pedal adjusts height from 36" to 62". Low profile feet slide easily under OR Table and hold stand in place.

Constructed of #304 stainless steel and all seams are fully welded. Brass inserts offer easy height adjustment and 2" **sayno2rust** composite casters provide ease of mobility. Ships fully assembled and ready to use.

- › Easy to clean
- › Sayno2rust composite casters
- › Low profile base

#### PROVEN DESIGN

Low profile feet slide easily under OR Table and hold stand in place.

#### MATERIAL

Constructed of #304 stainless steel and all seams are fully welded. Brass inserts offer easy height adjustment.

#### CASTERS

Are you tired of seeing rust on casters in your Operating Room? Our **sayno2rust** casters are standard on our Mayo Stands. They are non-marking neoprene twin wheel casters. In addition, these wheels are very easy to change with standard tools. They will make your OR run smoother and make your equipment last longer. **Improvinglife.**



Removable **tray** for easy cleaning



**Foot pedal** adjusts height from 36" to 62"

# Mayo Stands

## HCM760 Series

### Thumb Control Mayo Stands

Single post mayo stand with thumb control button. Removable tray for easy cleaning. Stainless steel thumb control adjusts height from 36" to 62". Low profile feet slide easily under OR Table and hold stand in place. Constructed of #304 stainless steel and all seams are fully welded. Brass inserts offer easy height adjustment and 2" **sayno2rust** composite casters provide ease of mobility. Ships fully assembled and ready to use.

- › **Can be controlled from anywhere without compromising the drape**
- › **No need to be standing to adjust stand**
- › **Low profile design slides under OR table base**

HCM760	12 5/8"W x 19 1/8"H Tray Size
HCM761	16 1/4"W x 21 1/4"H Tray Size

#### MOBILE BASE

Mobile Base can be added to any Mayo stand. Just add **"-MB"** to the SKU#.



Removable **tray** for easy cleaning



**Thumb control** adjusts height from 36" to 62"

# Mayo Stands

## HCM730 Series

### Knob Control Mayo Stands

Fully stainless steel welded construction but more economical version. Knob control with no-lose thumb knob. Designed for low volume OR's or physician's office. Removable tray for easy cleaning. Stainless steel thumb knob adjusts height from 36" to 62". Low profile feet slide easily under OR Table or exam table to hold stand in place. Constructed of #304 stainless steel and all seams are fully welded. Brass inserts offer easy height adjustment and 2" **sayno2rust** composite casters provide ease of mobility. Ships fully assembled and ready to use.

- › More economical but still a SS Mayo Stand
- › No-lose thumb knob
- › Low profile design to slide under OR Table or exam table

#### PROVEN DESIGN

Low profile feet slide easily under OR Table and hold stand in place.



#### MATERIAL

Constructed of #304 stainless steel and all seams are fully welded. Brass inserts offer easy height adjustment.

#### CASTERS

Are you tired of seeing rust on casters in your Operating Room? Our **sayno2rust** casters are standard on our Mayo Stands. They are non-marking neoprene twin wheel casters. In addition, these wheels are very easy to change with standard tools. They will make your OR run smoother and make your equipment last longer. **Improvinglife.**

# sayno2rust.<sup>®</sup>



Removable **tray** for easy cleaning



**Thumb knob** adjusts height from 36" to 62"

# Mayo Stands

## HCM770 Series

### Stainless Steel Ventric Stand

Built because of the request for a larger mayo stand, our lift assist Mayo Stand is the daddy of all Mayo Stands. Featuring a larger top, Crank lift, and large base, this is a heavy duty mayo stand. Top is 16 gauge material and is extremely rigid so when you put heavier objects on it, it doesn't deflect. Height adjusts from 36" to 56". Has two locking casters to hold it in place during surgery but easy to transport when needed. All edges including the handle are smooth to prevent puncture of the drape. Fully welded stainless steel construction. Ships fully assembled.

- › Ideal for large cases, especially ortho
- › Easy crank adjust to ergonomic height
- › Extra large surface and the needed stability

HCM770	30"W x 26"H Top Size
HCM771	32"W x 30"H Top Size



Large 4" **casters** with top brake for securing stand in place.



Extra large **top size** is over 50% larger than our largest mayo stand. Great for large cases, especially ortho



**Rail** around top to prevent instruments from rolling off top during surgery.



# Mayo Stands



## HCM700 Series Aluminum Base Mayo Stands

The Aluminum base mayo stand is great for physicians offices, emergency rooms or minor procedure rooms. Adjusts in height by a thumb knob below the top. Pole and top are stainless steel with aluminum base and twin rolling plastic casters.

- › Economical design yet durable
- › Removable Tray for easy cleaning
- › Two tray sizes available
- › Solid stainless steel pole construction
- › 2" twin wheel plastic casters
- › Knob operated lock adjusts height from 36"-62"

HCM705	12 5/8"W x 19 1/8"H Tray Size
HCM706	16 1/4"W x 21 1/4"H Tray Size



## HCM710 Series Space Saving Base Mayo Stands

The weighted base mayo stand is great for physicians offices, emergency rooms or minor procedure rooms. Adjusts in height by a thumb knob below the top. Pole and top are stainless steel. Epoxy coated base available in 5 colors to help each department maintain their assets. Simply add a - to the part number and follow with the first letter of the color..

- › Economical design yet durable
- › Removable Tray for easy cleaning
- › Two tray sizes available
- › Solid stainless steel pole construction
- › 3" neoprene non-marking casters
- › Knob operated lock adjusts height from 36"-62"

HCM710	12 5/8"W x 19 1/8"H Tray Size
HCM711	16 1/4"W x 21 1/4"H Tray Size

### Standard base colors are:

<span style="background-color: #333; width: 15px; height: 15px; display: block;"></span>	CHARCOAL	<span style="background-color: #C00; width: 15px; height: 15px; display: block;"></span>	R: RED	<span style="background-color: #000; width: 15px; height: 15px; display: block;"></span>	B: BLUE
<span style="background-color: #FFD700; width: 15px; height: 15px; display: block;"></span>	Y: YELLOW	<span style="background-color: #006400; width: 15px; height: 15px; display: block;"></span>	G: GREEN		



## HCM740 Series Mayo Tray Attachment

Mayo tray attaches directly to OR table rail for convenience and to save space in the OR. 5/8"mounting post attaches to the OR table rail with standard clark socket (*sold separately*). Height can be adjusted up to 26" above the table surface. 12 5/8" x 19 1/8" tray size.

HCM740	Mayo Tray Attachment
--------	----------------------

Attaches to rail with Clark Socket  
(*sold separately*)

**SEE PAGE 34**

# Additional OR Equipment



## HCM1000 Series Single Bowl Solution Stands

Stainless steel solution stand with one bowl. 16 gauge #304 stainless steel construction. Rubber bumpers inside basin ring keep basin firmly in place and reduce noise when rolling. 3" scuff resistant neoprene casters for ease of rolling. Measures 15"W x 18"L x 34"H. Available with shelf or H-brace depending on your preference. Ships fully assembled. 7 Qt. basin is included.

› Caster upgrade available

**sayno2rust.**

HCM1000	Single Bowl Solution Stand with H-Brace
HCM1001	Single Bowl Solution Stand with Shelf
HCM1007	Basin Only



## HCM1000 Series Double Bowl Solution Stands

Stainless Steel solution stand with two bowls. 16 gauge #304 stainless steel construction. Rubber bumpers inside basin ring keep basin firmly in place and reduce noise when rolling. 3" scuff resistant neoprene casters for ease of rolling. Measures 15"W x 33"L x 34"H. Available with shelf or H-brace depending on your preference. Ships fully assembled. 7 Qt. basin is included.

› Caster upgrade available

**sayno2rust.**

HCM1002	Double Bowl Solution Stand with H-Brace
HCM1003	Double Bowl Solution Stand with Shelf
HCM1007	Basin Only



## HCM540 Series Kick Bucket with Stand

Stainless steel Kick bucket with fully welded stainless steel frame. 2" swivel neoprene casters make it easy to move. Rubber bumpers hold bucket firmly in place. X-shaped base support keeps bucket supported off the floor even when it is full. Bucket has 12 qt. capacity and is removable. Perimeter rubber bumper protects walls. Ships fully assembled.

HCM540	Kick Bucket with Stand
HCM545	12 qt. Bucket Only



## HCM540 Series Sponge Receptacle with Stand

Stainless steel sponge receptacle with fully welded stainless steel frame. 2" swivel neoprene casters make it easy to move. Rubber bumpers hold basin firmly in place. X-shaped base support keeps basin supported off the floor even when it is full. Basin has 8.5 qt. capacity and is removable. Perimeter rubber bumper protects walls. Ships fully assembled.

HCM541	Sponge Receptacle with Stand
HCM542	8.5 qt. Receptacle only

# Additional OR Equipment



## HCM130 Series

### Chrome Step Stools

HCM138

Chrome Step Stool

Economic Step Stool for use anywhere you need an extra few inches. Very inexpensive when you compare to strains from over reaching and/or falls from standing on incorrect step stools such as chairs or counters. Constructed of chrome plated steel. Top platform measures 11 3/8" x 14 1/4" with anti-slip rubber tread. 9" platform height.



## HCM140 Series

### Stacking Interlocking Step Stool

HCM140

Stacking  
Interlocking  
Step Stool

Stacking interlocking step stools provide a great solution to varying heights of staff in the OR. These stools can be used individually or together to create a platform. They also can be used stacked on top of each other to create higher platforms. They are constructed of 16 gauge #304 stainless steel with the same operating room specs as all our other products. Slip proof rubber top for safety. Nest together easily with integral lip to hold them in place. Top size is 12" x 18" and height is 5".

- › Stack for easy storage
- › Great way to get everyone at the right height in the OR



## HCM140 Series

### Step Stool Carry Cart

HCM145

Step Stool  
Carry Cart

Carry Cart designed to hold 6 step stools. Ideal for storage and transport of step stools from one OR to another. Two locking casters secure cart in place. Fully welded stainless steel construction for easy cleaning.

- › Easy to clean
- › Saves time in transport and stores as well



## HCM140 Series

### Stirrup Carry Cart

HCM148

Stirrup Carry Cart

Cart designed specifically for the storage and transport of stirrups. Accessory rail to hold clamps and hooks hold most types of stirrups. Other accessories can be stored in lower tray. Constructed of 16 gauge stainless steel for easy cleaning.

- › Makes it easy to store and transport stirrups
- › Keeps clamps and stirrups together and clean

# Storage, Transfer and Transport

Intelligent solutions for the safe and efficient storage, transfer and transport of patients and supplies in healthcare facilities



# Storage, Transfer and Transport



## HCM2000 Series Round Stainless Steel Hampers

Constructed of 16 gauge fully welded stainless steel, these hampers will not rust. Lower shelf supports linen bag as it fills to prevent the bag dragging on the floor. Available in 2 different sizes to meet your needs. Fit standard size linen bags. 3" neoprene casters for ease of mobility even when full. 35" height for ease of loading.

› **Caster upgrade available**

**sayno2rust\***

HCM2000	18" Round Stainless Steel Hamper
HCM2001	25" Round Stainless Steel Hamper



## HCM2000 Series Triangle Stainless Steel Hampers

Constructed of 16 gauge fully welded stainless steel, these hampers will not rust. Lower shelf supports linen bag as it fills to prevent the bag dragging on the floor. Available in 2 different sizes to meet your needs. Fit standard size linen bags. 3" neoprene casters for ease of mobility even when full. 35" height for ease of loading.

› **Caster upgrade available**

**sayno2rust\***

HCM2006	18" Triangle Stainless Steel Hamper
HCM2005	25" Triangle Stainless Steel Hamper



## HCM2010 Series Square Chrome Hamper

Economical chrome hamper with flip top lid. Priced so you can have anywhere you need it. Fits 18" hamper bag. Plastic top can be easily cleaned. Foot control open for hands free operation. Hamper measures 19"W x 19 1/2"D x 37 3/4"H.

HCM2010	Square Chrome Hamper with Lid
---------	-------------------------------

# Storage, Transfer and Transport

## HCM3000 Series

### Utility Carts

These stainless steel utility carts are constructed of 16 gauge stainless steel. Fully welded, they are easily customized to meet your needs. Available standard with 2 or 3 shelves. Push handle is constructed of stainless steel. Entire cart is stainless steel for easy cleaning. Available with locking casters to hold cart in place when not in transit. 3" neoprene casters are standard. Ships fully assembled.

› Caster upgrade available

**sayno2rust.**

HCM3008	16"W x 24"L x 34"H, Two shelves
HCM3009	16"W x 24"L x 34"H, Three shelves
HCM3006	24"W x 30"L x 34"H, Two shelves
HCM3007	24"W x 30"L x 34"H, Three shelves
HCM3010	18"W X 27"L X 34"H, Two shelves
HCM3011	18"W x 27"L x 34"H, Three shelves

*Customize your cart:*

Choose your shelf size

Choose your height

Choose the qty of shelves  
and spacing



Additional  
shelves.



3" neoprene casters,  
locking casters  
available.



Stainless steel  
push handle cleans  
easily to minimize  
infection transfer.

# Storage, Transfer and Transport



## HCM120 Series

### Slider Style Patient Transfer Boards

Constructed of HDPE. 3/16" thick. Non-anti-static. Generous handle for easy transfers. Easily transfer patients from bed to gurney or OR Table. Smooth easy to clean surface. Available in two standard sizes.

HCM121	22"W x 33)L x 3/16"H
HCM122	22"W x 72)L x 3/16"H



## HCM120 Series

### Patient Transfer Roller Boards

Roller style patient transfer boards easily transfer patients from bed to gurney to OR table. Constructed of aluminum for easy to clean and minimal weight reasons. Standard style patient transfer board. 15" width and 3 lengths to meet your needs.

HCM123	15"W x 25)L
HCM125	15"W x 30)L
HCM130	15"W x 67)L

Replacement covers fit corresponding length of roller boards. Not guaranteed to fit other brands of roller boards.

HCM124	25)L Replacement Cover
HCM126	30)L Replacement Cover
HCM131	67)L Replacement Cover



## HCM190 Series

### EZ-Tranz Gold Patient Transfer Board

Gold Antimicrobial cover around a soft durable core, make this EZ-Tranz Board an elite transfer aid. Conveyor type action allows for smooth lateral transfer of patients. Extremely lightweight and easy to use make this the first choice for medical professionals and patients.

- › Easy to clean.
- › Antimicrobial
- › Extremely lightweight

Replacement covers fit corresponding length of transfer boards. Not guaranteed to fit other brands of transfer boards.

HCM190	15"W x 30)L
HCM191	15"W x 48)L

HCM195	30)L Replacement Cover
HCM196	48)L Replacement Cover

# Storage, Transfer and Transport



## HCM120 Series

### Wall Rack for Slider Style Patient Transfer Boards

HCM120-S

Wall Rack for Slider Style Patient Transfer Boards

Constructed of stainless steel material for ease of cleaning and prevention of infection spread. Hook holds small or large size patient transfer slider boards. Easy to hand and keeps boards clean and out of the way, yet in the spot you need them. Each rack holds one board. Mounting hardware not included due to varying substrates.



## HCM120 Series

### Wall Rack for Patient Transfer Roller Boards

HCM132

Wall Rack for Patient Transfer Roller Boards

Constructed of stainless steel material for ease of cleaning and prevention of infection spread. Fits all types of roller boards. Easy to hand and keeps boards clean and out of the way, yet in the spot you need them. Each rack holds one board. Mounting hardware not included due to varying substrates.



## HCM190 Series

### Wall Rack for EZ-Tranz Gold Patient Transfer Boards

HCM199

Wall Rack for EZ-Tranz Gold Patient Transfer Boards

Constructed of stainless steel material for ease of cleaning and prevention of infection spread. Fits all sizes of EZ-Tranz Gold transfer boards. Keeps boards clean and out of the way, yet in the spot you need them. Each rack holds one board. Mounting hardware not included due to varying substrates.

# Storage, Transfer and Transport



## HCM670 Series

### Glove Boxes

Constructed of durable stainless steel. Easy to clean. Manufactured to our OR specifications for minimizing glove tears and not harboring bacteria. Durable one-piece construction. Available in horizontal and vertical options.

HCM670	Single Box
HCM671	Double Box Vertical
HCM672	Double Box Horizontal
HCM673	Triple Box Vertical
HCM674	Triple Box Horizontal
HCM675	Quad Box Vertical



## HCM640 Series

### Wall Shelves

Stainless steel 16 gauge one piece construction. Constructed to our OR specs for minimizing of bacteria and easy to clean. Available in 8", 10", 12", 15" and 18" deep. All shelves over 5' included center support. Mounting hardware not included due to differences in substrate. Custom sizes available.

HCM640	24" Wide x 10" Deep Stainless Steel Wall Shelf	HCM664	72" Wide x 12" Deep Stainless Steel Wall Shelf
HCM641	30" Wide x 10" Deep Stainless Steel Wall Shelf	HCM665	84" Wide x 12" Deep Stainless Steel Wall Shelf
HCM642	36" Wide x 10" Deep Stainless Steel Wall Shelf	HCM666	96" Wide x 12" Deep Stainless Steel Wall Shelf
HCM643	48" Wide x 10" Deep Stainless Steel Wall Shelf	HCM667	30" Wide x 12" Deep Stainless Steel Wall Shelf
HCM644	60" Wide x 10" Deep Stainless Steel Wall Shelf	HCM680	24" Wide x 15" Deep Stainless Steel Wall Shelf
HCM645	72" Wide x 10" Deep Stainless Steel Wall Shelf	HCM681	30" Wide x 15" Deep Stainless Steel Wall Shelf
HCM646	84" Wide x 10" Deep Stainless Steel Wall Shelf	HCM682	36" Wide x 15" Deep Stainless Steel Wall Shelf
HCM647	96" Wide x 10" Deep Stainless Steel Wall Shelf	HCM683	48" Wide x 15" Deep Stainless Steel Wall Shelf
HCM650	24" Wide x 8" Deep Stainless Steel Wall Shelf	HCM684	60" Wide x 15" Deep Stainless Steel Wall Shelf
HCM651	36" Wide x 8" Deep Stainless Steel Wall Shelf	HCM685	72" Wide x 15" Deep Stainless Steel Wall Shelf
HCM652	48" Wide x 8" Deep Stainless Steel Wall Shelf	HCM686	84" Wide x 15" Deep Stainless Steel Wall Shelf
HCM653	60" Wide x 8" Deep Stainless Steel Wall Shelf	HCM687	96" Wide x 15" Deep Stainless Steel Wall Shelf
HCM654	72" Wide x 8" Deep Stainless Steel Wall Shelf	HCM690	24" Wide x 18" Deep Stainless Steel Wall Shelf
HCM655	84" Wide x 8" Deep Stainless Steel Wall Shelf	HCM691	30" Wide x 18" Deep Stainless Steel Wall Shelf
HCM656	96" Wide x 8" Deep Stainless Steel Wall Shelf	HCM692	36" Wide x 18" Deep Stainless Steel Wall Shelf
HCM657	30" Wide x 8" Deep Stainless Steel Wall Shelf	HCM693	48" Wide x 18" Deep Stainless Steel Wall Shelf
HCM660	24" Wide x 12" Deep Stainless Steel Wall Shelf	HCM694	60" Wide x 18" Deep Stainless Steel Wall Shelf
HCM661	36" Wide x 12" Deep Stainless Steel Wall Shelf	HCM695	72" Wide x 18" Deep Stainless Steel Wall Shelf
HCM662	48" Wide x 12" Deep Stainless Steel Wall Shelf	HCM696	84" Wide x 18" Deep Stainless Steel Wall Shelf
HCM663	60" Wide x 12" Deep Stainless Steel Wall Shelf	HCM697	96" Wide x 18" Deep Stainless Steel Wall Shelf

# Accessories for OR Tables

Accessories that attach to the OR Table to support and hold patients' anatomy during surgery for optimum surgical results.



# Accessories for OR Tables



Underside View

## HCM401 Series

### Radiolucent Height Adjustable Armboard

This armboard is designed to support patient's arms during surgery. Large portion of board is also radiolucent for C-arm capabilities. Rotates 180 degrees to desired position. It also adjusts in height to match the top of the OR table or to elevate the patient's arm if needed. Adjusts up to 4" in height. Standard size is 6" x 26". Clamps securely to any American made OR table rail.

- › Radiolucent for C-arming
- › Adjusts in height to accommodate thicker OR table pads

HCM401

Radiolucent  
Armboard

Uses 26" L Armboard Pad  
Pads sold separately

[SEE PAGE 33](#)



Underside View

## HCM410 Series

### Droplatch Armboard with 2" Pad

This armboard is designed to support patient's arms during surgery. Standard droplatch style armboard for ease of use. Rotates 180 degrees to desired position with lever lock. Available in either 24" or 26" length. Clamps to any American made OR table rail. Includes a 2" pad.

- › Easy to use and adjust
- › Available in two lengths

HCM410-  
W2P

24" L Drop-Latch  
Armboard with 2" Pad

HCM411-  
W2P

26" L Drop-latch  
Armboard with 2" Pad



Underside View

## HCM410 Series

### Droplatch Armboard without Pad

This armboard is designed to support patient's arms during surgery. Standard droplatch style armboard for ease of use. Rotates 180 degrees to desired position with lever lock. Available in either 24" or 26" length. Clamps to any American made OR table rail.

- › Easy to use and adjust
- › Available in two lengths

HCM410

24" L Drop-Latch  
Armboard

HCM411

26" L Drop-latch  
Armboard

*Pads sold separately*

[SEE PAGE 33](#)

# Accessories for OR Tables




HCM412 with 2" pad

## HCM410 Series Bariatric Droplatch Armboards

This armboard is designed to support bariatric patient's arms during surgery. Standard droplatch style armboard for ease of use. Extra wide to accommodate larger arms. Rotates 180 degrees to desired position with lever lock. Clamps to any American made OR table rail.

- › **Designed for Bariatric patients**
- › **10" width**

HCM412	26" L x 10" W Drop-Latch Bariatric Armboard
HCM413	26" L x 10" W x 2" H Bariatric Armboard Pad

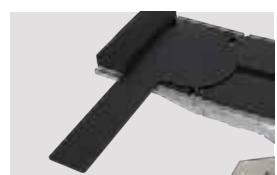

HCM440 under pad

## HCM440 Series Underpad Armboard

This armboard is designed to slide under the OR Table pad and be held in place by the weight of the patient. Can be adjusted as needed. Phenolic material is completely radiolucent for C-arm capabilities. Exposed Armboard portion is 6" W x 24" L.

- › **Radiolucent**
- › **Can be slid under the pad of any OR table or stretcher**

HCM440	Under Pad Armboard
--------	-----------------------

*Pads sold separately*

**SEE PAGE 33**



HCM442 with 2" pad

## HCM440 Series Extra Wide Underpad Armboard

This armboard is designed to slide under the OR Table pad and be held in place by the weight of the patient. Can be adjusted as needed. Phenolic material is completely radiolucent for C-arm capabilities. Exposed Armboard portion is 13.25" W x 24" L.

- › **Radiolucent**
- › **Can be slid under the pad of any OR table or stretcher**

HCM442	Extra Wide Under Pad Armboard
HCM443	21" L x 13.25" W x 2" H Extra Wide Armboard Pad

# Accessories for OR Tables



Underside View

## HCM410 Series

### Rail Locking Drop-Latch Style Armboard

This tried and proven armboard design is a great fit wherever you need an armboard. Easily attaches to any US made operating table side rail

- **Armboard drops onto rail and locks securely to the rail with a knob.**
- **Finger latch trigger on outer end of board allowing for 180 degree rotation**
- **Durable steel casting for longevity**
- **Armboard is made of 3/8" thick black phenolic and weighs a slim 4.5 pounds**
- **Pad sold separately**
- **Available in 24" or 26" long sizes**

HCM410-RL	24" L Rail-Lock Armboard
HCM411-RL	26" L Rail-Lock Armboard

*Pads sold separately*

**SEE PAGE 33**



HCM405 with pad


## HCM400 Series

### Width Extender and Pad

The OR table Width Extender is designed to add width to your OR table for those patients that are a little wider, or for specific procedures. Integral posi-lock rail clamps and aluminum side rail make this unit all-in-one. Available in 4" or 8" widths.

- **Equipped with integral "posi-lock" rail clamps and aluminum side rail**
- **Available in 4" or 8" widths**
- **Lock securely to existing siderail**
- **Quick set-up**
- **Pad sold separately**

HCM405	OR Table Width Extender 4" W x 20" H
HCM404	2" W Pad for Width Extender (HCM405)
HCM406	Bariatric OR Table Width Extender 8" W x 20" H
HCM407	2" W Pad for Bariatric Width Extender (HCM406)

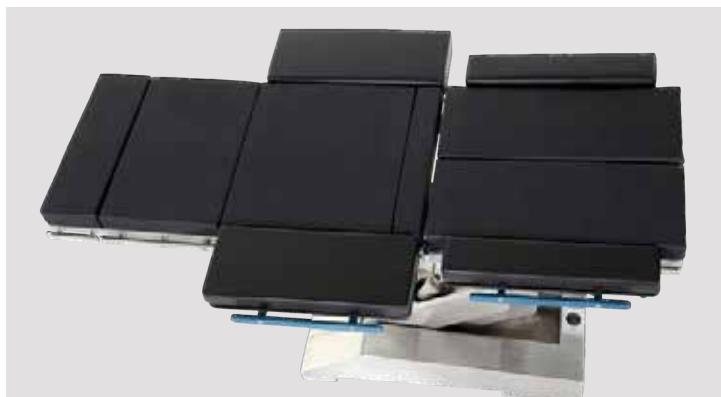

Add the OR Width Extender in the best configuration to suit the requirements of each individual procedure.

Shown here is an OR Table with HCM406 on either side of the torso and HCM405 on either side of the legs.

The possibilities are endless!

# Accessories for OR Tables



## HCM450 Series Elbow Arthroscopy Positioner

This surgical arm prepper places the patients elbow in optimal positioning for elbow arthroscopy. Easy access to a tourniquet without interfering with the sterile field. Prevents slipping without the use of straps. Durable Stainless steel construction Includes pressure diffusion pad to reduce impact on nerves and tissue.

HCM450	Elbow Arthroscopy Positioner with Pad
HCM451	Replacement Pad only for HCM450



## HCM440 Series Wall Mount Armboard Holders

Stainless steel rack holds armboards on the wall when not in use. Keeps them in a safe place to minimize damage and keep them where you need them. Constructed from 16 gauge stainless steel for easy cleaning and longevity. Will hold one armboard with standard 2" pad. Available for thicker pads by request.

- › Keeps armboards from getting dirty
- › Easy to clean

HCM445	Wall Mount Armboard Holder
HCM446	Wall Mount Bariatric Armboard Holder



## HCM400 Series Armboard Pads

Pads for armboards in your desired thickness and size. Standard foam with electrolyte cover like your OR table pads.

- › Easy to clean
- › Available in any thickness to match your OR table

HCM415	24" L x 2" H Armboard Pad
HCM416	26" L x 2" H Armboard Pad
HCM417	24" L x 3" H Armboard Pad
HCM418	26" L x 3" H Armboard Pad
HCM409	24" L x 4" H Armboard Pad
HCM419	26" L x 4" H Armboard Pad
HCM413	26" L x 10" W x 2" H Bariatric Armboard Pad
HCM443	21" L x 13.25" W x 2" H Extra Wide Armboard Pad

# Accessories for OR Tables



## HCM160 Series

### Stainless Steel Elite Clark Socket

HCM159

Stainless Steel  
Elite Clark Socket

Designed to clip on to the rail anywhere, not just at ends or at rail slots. Holds flat blade accessories as well as 3/4" round bar accessories with the same clamp. Solid stainless steel construction for strength and durability. Compatible with all American Made OR tables.

- › Two clamps in one
- › Clips onto rail anywhere

.....



## HCM160 Series

### Stainless Steel Clark Socket

HCM160

Stainless Steel  
Clark Socket

Designed to hold round post accessories up to 11/16" diameter. Solid stainless steel construction for strength and durability. Compatible with American made OR Table rails. Clips onto rail at notches or rail ends.



## HCM160 Series

### Stainless Steel Side Rail Socket

HCM162

Stainless Steel Side  
Rail Socket

Designed to hold flat blade accessories that measure 1" x 1/4". Solid stainless steel construction for strength and durability. Compatible with American made OR Table rails. Clips onto rail at rail ends.



## HCM160 Series

### Aluminum Side Rail Socket

HCM161

Aluminum Side  
Rail Socket

Designed to hold flat blade accessories that measure 1" x 1/4". Aluminum construction make this clamps lighter but yet durable. Compatible with American made OR Table rails. Clips onto rail at notches or rail ends.



## HCM160 Series

### Universal Rail Clamp

HCM163

Universal  
Rail Clamp

Designed to hold round and flat blade accessories. Aluminum construction makes this clamp lighter but yet durable. Compatible with American made OR Table rails. Clamps onto rail anywhere.

# Accessories for OR Tables



## HCM100 Series Anesthesia Shield

HCM107 | Stainless Steel Anesthesia Shield

Mounts to rail with two Side Rail Sockets (*sold separately*)

[SEE PAGE 34](#)



## HCM100 Series Rigid Anesthesia Screen

HCM105 | Rigid Anesthesia Screen

Mounts to rail with Clark Socket (*sold separately*)

[SEE PAGE 34](#)



## HCM100 Series Adjustable Anesthesia Screen

HCM101 | Adjustable Anesthesia Screen

Mounts to rail with Clark Socket (*sold separately*)

[SEE PAGE 34](#)



## HCM100 Series Malleable Anesthesia Screen

HCM106 | Malleable Anesthesia Screen

Mounts to rail with two Clark Sockets (*sold separately*)

[SEE PAGE 34](#)



## HCM100 Series Stainless Steel Hose Tree

HCM110 | Stainless Steel Hose Tree

Constructed of durable stainless steel, this hose tree is sure to last a long time. Easy to clean stainless steel. Efficiently keeps tubes and hoses in place and organized during surgery.

# Accessories for OR Tables



## HCM150 Series Stand Alone Leg Prep

The Stand Alone Leg Prep is adjusted completely independent from the OR Table and the base can slide in under the OR Table if needed. Perfect for cases when surgery is performed on a stretcher or gurney instead. Upper portion swivels 360 degrees for correct positioning. Adjusts in height with a thumb knob from 40"-60". Supports the legs of patients who are up to 300 lbs.

HCM154

Stand Alone Leg Prep



## HCM150 Series Single Loop Leg Prep with Pad

Stainless steel 5/8" diameter mounting post fits into standard clark socket (*sold separately*). Adjusts in height up to 30". Wipe-clean polyurethane pad. Weight limit is 50 lbs.

HCM157

Single Loop Leg Prep with Pad

HCM156

Replacement Pads Case of 12

Fits into standard Clark Socket (*sold separately*)

**SEE PAGE 34**



## HCM150 Series Candy Cane Stirrup

Constructed from 304 stainless steel for durability and longevity. Integral clamp allows height adjustment for 23"-41". Large 8" arc allows for additional leg clearance. Mounts to OR table rail with clark socket (*sold separately*). 1" wide cotton strap that conforms to 9" loops is included.

HCM150

Candy Cane Stirrup

HCM151

Replacement Cotton Straps

Mounts to rail with Clark Socket (*sold separately*)

**SEE PAGE 34**

# Accessories for OR Tables



## HCM600 Series Picket Fence Leg Holder

The Picket Fence Leg Prep Unit efficiently preps both legs at once. It has a 5/8" diameter, stainless steel support bar with an easily adjustable height up to 30". The Picket Fence Leg Prep Unit is easy to clean. Attaches to OR Table rail with clark socket (*sold separately*).

HCM611 Picket Fence Leg Holder

Mounts to rail with Clark Socket  
(*sold separately*)

**SEE PAGE 34**



## HCM600 Series Double Picket Fence Leg Holder

This double picket fence leg prepper efficiently preps both legs at once. Double flat blade supports bring added stability to this unit. Can adjust in height up to 30". Top is covered in rubber and lower supports are stainless steel making it easy to clean. Attaches to OR table rail with two side rail clamps (*sold separately*).

HCM610 Double Picket Fence Leg Holder

Attaches to two side rail clamps  
(*sold separately*)

**SEE PAGE 34**



## HCM610 Series TKR Support with Pad

Helps support legs for total knee replacement procedures. Constructed of stainless steel and polyoxy-methylene for ultimate support. Top portion is radiolucent. 5/8" diameter stainless steel post fits into clark socket (*sold separately*). Black pad included.

HCM612 TKR Support with Pad

HCM613 Replacement Pads Case of 12

Fits into Clark Socket  
(*sold separately*)

**SEE PAGE 34**

# Accessories for OR Tables



## HCM740 Series Mayo Tray Attachment

HCM740      Mayo Tray Attachment

Attaches to rail with Clark Socket  
(sold separately)

[SEE PAGE 34](#)



## HCM160 Series Insta-Rail

HCM168      Insta-Rail

Quickly and easily attach an OR table rail to any surface. An imaging table or an arm/hand table or any other surface up to 1 1/2" thick. Stainless steel rail is 6" long and standard American made rail dimensions. Holds any clamp that is compatible with American made OR table rails.



## HCM160 Series Side Rail Extension

HCM167      Side Rail Extension

Need to have a little longer side rail? Look no further than this side rail extension. Clamps on to the OR table rail and extends it 6". Compatible with any American made OR table rail.



## HCM140 Series Accessory and Stirrup Carry Cart

HCM148      Accessory and Stirrup Carry Cart

Cart designed specifically for the storage and transport of stirrups and other OR table accessories. Accessory rail to hold clamps and hooks hold most types of stirrups. Other accessories can be stored in lower tray. Constructed of 16 gauge stainless steel for easy cleaning.

- › **Makes it easy to store and transport stirrups**
- › **Keeps clamps and stirrups together and clean**

# Arm and Hand Surgery Tables

Surfaces that attach to the OR table or stretcher for performing arm and hand surgeries.



# Arm and Hand Surgery Tables

## Arm and Hand Tables without Leg Support

### HCM300 Series

#### Free Access Arm and Hand Table



This one-of-a-kind shape will make any surgeon smile. It allows easy access to their work with its narrow design and no leg to obstruct c-arm or X-ray procedures. Radiolucent center allows c-arm access. Clamps securely to the OR table rail with our integral posi-lock clamps. No extra clamps needed or to get misplaced. Our most compact arm/hand surgery table with 12" x 30" operating surface. Great for many procedures including fistula, angiogram and angioplasty.

HCM301	Free Access Arm and Hand Table with 2" Pad
HCM302	Replacement 2" Pad



### HCM330 Series

#### Minor Procedure Arm and Hand Table

Phenolic material for longevity and durability  
15" x 24" Minor Procedure size  
Integral Clamps don't get lost and make it easy to clamp to the OR Table side rail  
No obstruction to C-arm access  
2" pad included, other thicknesses available to suit your application.  
Facilitates arm or hand surgery on Any US Made OR Table

HCM333	Minor Procedure Arm and Hand Table with 2" Pad
HCM334	Replacement 2" Pad



### HCM330 Series

#### Carpal Arm and Hand Table

This rectangle shaped arm/hand table offers superior rigidity without the obstruction of a lower leg. Radiolucent center portion allows c-arm access. Clamps securely to the OR table rail with our integral posi-lock clamps. No extra clamps needed or to get misplaced. Full size 18" x 32" operating surface and plenty of lower support. This is the ultimate choice for carpal tunnel procedures and many others.

HCM330	Carpal Arm and Hand Table with 2" Pad
HCM312	Replacement 2" Pad



### HCM310 Series

#### Carbon Fiber Rectangle Arm and Hand Table

This rectangle shaped arm/hand table made of lightweight but durable carbon fiber. Completely radiolucent for excellent C-arm image. Clamps securely to the OR table rail with our integral posi-lock clamps. No extra clamps needed or to get misplaced. Full size 18" x 32" operating surface but yet only weighs a mere 9 lbs. Perfect choice for minor arm/hand procedures.

HCM318	Carbon Fiber Rectangle Arm and Hand Table with 2" Pad
HCM312	Replacement 2" Pad

# Arm and Hand Surgery Tables

## Arm and Hand Tables with Leg Support That Clamp to the OR Table Rail



### HCM310 Series

#### Rectangle Phenolic Arm and Hand Tables

This rectangle arm/hand table is constructed of durable phenolic like the surface of the OR table. It clamps securely to the OR table rail and has a drop down leg for added support. Radiolucent center portion allows c-arm access. Clamps securely to the OR table rail with our integral posi-lock clamps. No extra clamps needed or to get misplaced. Full size 18" x 32" operating surface and plenty of lower support. Leg adjusts up to 12" to meet the ideal height for the surgeon.

HCM310	Rectangle Phenolic Arm and Hand Table, Single Foot, with 2" Pad
HCM311	Rectangle Phenolic Arm and Hand Table, Double Foot, with 2" Pad
HCM312	Replacement 2" Pad



### HCM310 Series

#### Rectangle Carbon Fiber Arm and Hand Tables

Constructed of durable and lightweight carbon fiber, this table has superior c-arm image as well as being very lightweight, but is the strongest table available. It clamps securely to the OR table rail and has a drop down leg for added support. Radiolucent center portion allows c-arm access. Clamps securely to the OR table rail with our integral posi-lock clamps. No extra clamps needed or to get misplaced. Full size 18" x 32" operating surface and plenty of lower support. Leg adjusts up to 12" to meet the ideal height for the surgeon.

HCM315	Rectangle Carbon Fiber Arm and Hand Table, Single Foot, with 2" Pad
HCM316	Rectangle Carbon Fiber Arm and Hand Table, Double Foot, with 2" Pad
HCM312	Replacement 2" Pad



### HCM320 Series

#### Hourglass Phenolic Arm and Hand Tables

This hourglass shaped arm/hand table is constructed of durable phenolic like the surface of the OR table. The hourglass shape extends to 24" width to provide wrist support for the surgeon. It clamps securely to the OR table rail and has a drop down leg for added support. Radiolucent center portion allows c-arm access. Clamps securely to the OR table rail with our integral posi-lock clamps. No extra clamps needed or to get misplaced. 15" x 32" operating surface and plenty of lower support. Leg adjusts up to 12" to meet the ideal height for the surgeon.

HCM320	Hourglass Phenolic Arm and Hand Table, Single Foot, with 2" Pad
HCM321	Hourglass Phenolic Arm and Hand Table, Double Foot, with 2" Pad
HCM322	Replacement 2" Pad



### HCM320 Series

#### Hourglass Carbon Fiber Arm and Hand Tables

Constructed of durable and lightweight carbon fiber, this table has superior c-arm image as well as being very lightweight, but is the strongest table available. The hourglass shape extends to 24" width to provide wrist support for the surgeon. It clamps securely to the OR table rail and has a drop down leg for added support. Radiolucent center portion allows c-arm access. Clamps securely to the OR table rail with our integral posi-lock clamps. No extra clamps needed or to get misplaced. Full size 15" x 32" operating surface and plenty of lower support. Leg adjusts up to 12" to meet the ideal height for the surgeon.

HCM325	Hourglass Carbon Fiber Arm and Hand Table, Single Foot, with 2" Pad
HCM326	Hourglass Carbon Fiber Arm and Hand Table, Double Foot, with 2" Pad
HCM322	Replacement 2" Pad

# Arm and Hand Surgery Tables

## Arm and Hand Tables with Outer Leg That Slide Under The Pad

As Ambulatory Surgery Centers become more common in the medical field, many surgeons wish to perform Arm/Hand surgery on a Stretcher or Stretcher Chair. We have an extended range of Arm/Hand tables that slide under the pad of these units and are held in place by the weight of the patient. All tables include 2" thick standard pad.



### HCM330 Series

#### Rectangle Phenolic Under Pad Arm and Hand Tables

Our most popular size and shape of Arm/Hand table is the rectangle. Large 18" x 32" operating surface with excellent C-arm capabilities make this a first choice. Phenolic is our most popular material for arm/hand surgery tables. Durable with good C-arm capabilities it makes a great table. Available with drop down leg and either single or double foot.

HCM335	Phenolic, Rectangle Under Pad Arm and Hand Table with Single Foot
HCM336	Phenolic, Rectangle Under Pad Arm and Hand Table with Double Foot

HCM336  
on table



HCM336  
underside  
folded



### HCM330 Series

#### Hourglass Phenolic Under Pad Arm and Hand Tables

This hourglass shaped arm/hand table is constructed of durable phenolic like the surface of the OR table. The hourglass shape extends to 24" width to provide wrist support for the surgeon. It slides easily under the pad and is held in place by the patients weight. Radiolucent center portion allows c-arm access. 15" x 32" operating surface and plenty of lower support. Leg adjusts up to 12" to meet the ideal height for the surgeon.

HCM337	Phenolic, Hourglass Under Pad Arm and Hand Table with Single Foot
HCM338	Phenolic, Hourglass Under Pad Arm and Hand Table with Double Foot

HCM338  
on table



HCM338  
underside  
folded



# Arm and Hand Surgery Tables



Single-Foot

HCM340  
on table



HCM340  
underside  
folded



HCM340

Stainless Steel,  
Rectangle Under Pad  
Arm and Hand Table  
with Single Foot

HCM341

Stainless Steel,  
Rectangle Under Pad  
Arm and Hand Table  
with Double Foot



Single-Foot

HCM343  
on table



HCM343  
underside  
folded

HCM342

Carbon Fiber,  
Rectangle Under Pad  
Arm and Hand Table  
with Single Foot

HCM343

Carbon Fiber,  
Rectangle Under Pad  
Arm and Hand Table  
with Double Foot

# Arm and Hand Surgery Tables

## Arm and Hand Tables with Mobile Base

Arm hand tables with Mobile bases offer the best of versatility. Able to roll from room to room without the nurses having to lift them is a great benefit. Adjusts in height from 30"-42". They come equipped with two locking casters and two swivel casters. Several different shape and size options of tops ensure your needs are met.



HCM310-MB  
Underside View

### HCM310 Series

#### Rectangle Arm and Hand Tables with Mobile Base

Our most popular size and shape of Arm/Hand table is the rectangle. Large 18" x 32" operating surface with excellent c-arm capabilities make this a first choice. Several different options are available to secure this unit to the table or stretcher.

**Roll Up:** This option simply rolls up to the OR table or stretcher and is held in place by the locking casters. Works for performing surgery on OR table, stretcher and/or stretcher chair.

**Clamps:** Includes our integral posi-lock clamps to ensure that the table remains securely attached to the OR table during surgery. A great option if you are performing surgery strictly on an OR table.

**Underpad:** Has a longer portion that can slide under the pad of the OR table or stretcher and is held in place by the patients weight. Works for performing surgery on OR table, stretcher and/or stretcher chair.

HCM310-MB	Rectangle Arm and Hand Table with Mobile Base, 2" Pad Included
HCM310-MBCL	Rectangle Arm and Hand Table with Mobile Base, Clamps Included, 2" Pad Included
HCM310-MBUP	Rectangle Arm and Hand Table with Mobile Base, Slides Under Pad, 2" Pad Included



HCM310-MB  
Underside View

### HCM320 Series

#### Hourglass Arm and Hand Tables with Mobile Base

Hourglass shape is popular because it allows a wide 24" outer edge for wrist support of the surgeon during surgery. 15" x 32" operating surface with excellent c-arm capabilities make this table a great fit where an hourglass shape is preferred. Several different options are available to secure this unit to the table or stretcher.

**Roll Up:** This option simply rolls up to the OR table or stretcher and is held in place by the locking casters. Works for performing surgery on OR table, stretcher and/or stretcher chair.

**Clamps:** Includes our integral posi-lock clamps to ensure that the table remains securely attached to the OR table during surgery. A great option if you are performing surgery strictly on an OR table.

**Underpad:** Has a longer portion that can slide under the pad of the OR table or stretcher and is held in place by the patients weight. Works for performing surgery on OR table, stretcher and/or stretcher chair.

HCM320-MB	Hourglass Arm and Hand Table with Mobile Base with 2" Pad Included
HCM320-MBCL	Hourglass Arm and Hand Table with Mobile Base, Clamps Included, 2" Pad Included
HCM320-MBUP	Hourglass Arm and Hand Table with Mobile Base, Slides Under Pad, 2" Pad Included

# Arm and Hand Surgery Tables


## Arm and Hand Table Accessories



### HCM350 Series Horizontal Traction Towers

This traction tower bolts on to existing holes on end of Arm/Hand table brackets. Manufactured of quality stainless steel construction. 3/4" outside dimension tube with 1/2" inner tube that telescopes to desired position. All necessary hardware included.

HCM344	18" Horizontal Traction Tower for Rectangle Tables
HCM345	24" Horizontal Traction Tower for Hourglass Tables

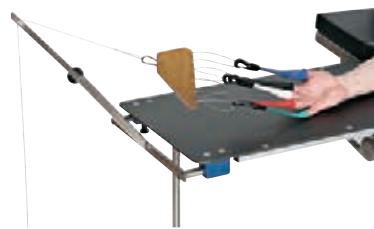

### HCM350 Series Equalizing Device with Finger Traps

Used to keep pressure on digits when 3 or more are in traction. Can be used with IV Poles or other extremity systems. Stainless steel hooks and chain are fully autoclavable. Five different sizes of finger traps available to meet your needs.

HCM350	Equalizing Device and Finger Traps
HCM355	X-small Finger Trap
HCM356	Small Finger Trap
HCM357	Medium Finger Trap
HCM358	Large Finger Trap
HCM359	X-large Finger Trap



### HCM350 Series Add-A-Rail

Need to add a rail onto the end of your Arm/Hand table for accessories to attach to? Solid stainless steel construction. Compatible to all American made OR table rail sockets. Easily mounts to end of Arm/Hand Table.

HCM314	18" Long Add-A-Rail for Rectangle Tables
HCM324	24" Long Add-A-Rail for Hourglass Table



### HCM300 Series Replacement Pads

Standard 2" pads are available for all our Arm/Hand Tables. Custom thicknesses are available in order to meet the thickness of your current OR table pads.

HCM302	2" Thick Standard Pad for Free Access Table
HCM312	2" Thick Standard Pad for Rectangle Table
HCM322	2" Thick Standard Pad for Hourglass Table

# Arm and Hand Surgery Tables

Item #	Size			Material			Base				Table Attachment		Pad Included
	Rectangle	Hourglass	Free Access	Carbon Fiber	Stainless Steel	Phenolic	Mobile Base	Flip Up Double Leg	Flip Up Single Leg	No Leg	OR Table Rail	Under Pad	
HCM301			X			X				X	X		X
HCM310	X					X			X		X		X
HCM310-MB	X					X	X						X
HCM310-MBCL	X					X	X				X		X
HCM310-MBUP	X					X	X					X	X
HCM311	X					X		X			X		X
HCM315	X			X					X		X		X
HCM316	X			X				X			X		X
HCM318	X			X						X	X		X
HCM320		X				X			X		X		X
HCM320-MB		X				X	X						X
HCM320-MBCL		X				X	X				X		X
HCM320-MBUP		X				X	X					X	X
HCM321		X				X		X			X		X
HCM325		X		X					X		X		X
HCM326		X		X				X			X		X
HCM330	X					X				X	X		X
HCM333	X					X				X	X		X
HCM335	X					X			X			X	X
HCM336	X					X		X				X	X
HCM337		X				X			X			X	X
HCM338		X				X		X				X	X
HCM340	X				X				X			X	X
HCM341	X				X			X				X	X
HCM342	X			X					X			X	X
HCM343	X			X				X				X	X



## IV Poles and Accessories

Slender Stands and accessories that assist in the suspension of containers of liquid for the administration of IV to patients.

# IV Poles and Accessories



## HCM210 Series

### 5 Leg Economy IV Poles

This 5 Leg IV Pole was designed with functionality in mind. The 5 leg design and 24" diameter base brings the needed stability. Neoprene 2" casters offer adequate mobility and are silent while transporting. 1" diameter pole is available in chrome finish or durable stainless steel finish. Pole adjusts in height from 52"-93" with no-lose friction knob. Rubber inserts on outer edges of legs protect walls and other objects from scuffing.


- **Various top options are available**

#### Chrome Poles

HCM201	Chrome IV Pole W/Thumb Knob, 2 Hook Top, 5 Leg Base W/2" Casters
HCM202	Chrome IV Pole W/Thumb Knob, 4 Hook Top, 5 Leg Base W/2" Casters

#### Stainless Steel Poles

HCM210	SS IV Pole W/2 Hook Top, No Lose Thumb Knob, 5 Leg Base W/2" Casters
HCM210-RT	SS IV Pole W/8 Hook Rake-Top, No Lose Thumb Knob, 5 Leg Base W/2" Casters
HCM211	SS IV Pole W/4 Hook Top, No Lose Thumb Knob, 5 Leg Base W/2" Casters
HCM212	SS IV Pole W/6 Hook Top, No Lose Thumb Knob, 5 Leg Base W/2" Casters

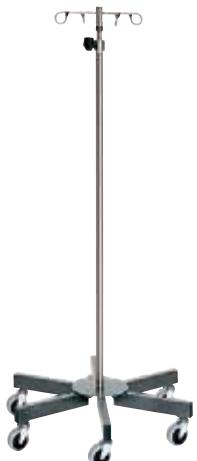

## HCM230 Series

### 6 Leg Economy IV Poles

This 6 Leg IV Pole was designed with functionality in mind. The 6 Leg design and 24" diameter base brings the needed stability. Neoprene 3" casters offer ease of mobility and are silent while transporting. 1" diameter pole is available in chrome finish or durable stainless steel finish. Pole adjusts in height from 52"-93" with no-lose friction knob. Rubber inserts on outer edges of legs protect walls and other objects from scuffing.

- **Various top options are available**
- **Base is epoxy coated in the standard charcoal color**

#### Chrome Poles

HCM234	Chrome IV Pole W/Thumb Knob, 2 Hook Top, 6 Leg Base W/3" Casters
HCM235	Chrome IV Pole W/Thumb Knob, 4 Hook Top, 6 Leg Base W/3" Casters

#### Stainless Steel Poles

HCM237	SS IV Pole W/2 Hook Top, No Lose Thumb Knob, 6 Leg Base W/3" Casters
HCM237-RT	SS IV Pole W/8 Hook Rake-Top, No Lose Thumb Knob, 6 Leg Base W/3" Casters
HCM238	SS IV Pole W/4 Hook Top, No Lose Thumb Knob, 6 Leg Base W/3 Casters
HCM239	SS IV Pole W/6 Hook Top, No Lose Thumb Knob, 6 Leg Base W/3" Casters



## HCM220 Series

### 5 Leg Space-Saving IV Poles

This 5 Leg Space-saving IV Pole was designed to save valuable space in the hospital without sacrificing stability. Base is 16" diameter and weighs 22 lbs. 30 scuff resistant casters are standard for ease of mobility. Height adjusts from 52"-94" with no-lose thumb knob. Available in chrome finish or durable stainless steel finish. Different top options are available to meet your needs. Designed to hold 1-2 infusion pumps. Standard epoxy coated base color is charcoal.

Colored bases are available at no extra charge to help each department maintain their assets. Simply add a ~ to the part number and follow with the first letter of the color.

Standard colors are:

CHARCOAL	R: RED	B: BLUE
Y: YELLOW	G: GREEN	

#### Chrome Poles

HCM220-2	Chrome IV Pole W/Thumb Knob, 2 Hook Top, 5 Leg 22 lb. Base W/3" Casters, Charcoal
HCM220	Chrome IV Pole W/Thumb Knob, 4 Hook Top, 5 Leg 22 lb. Base W/3" Casters, Charcoal

#### Stainless Steel Poles

HCM221	SS IV Pole W/2 Hook Top, No Lose Thumb Knob, 5 Leg 16"dia. 22 lb. Base W/3" Casters, Charcoal
HCM222-RT	SS IV Pole W/8 Hook Rake-Top, No Lose Thumb Knob, 5 Leg 16"dia. 22 lb. Base W/3" Casters, Charcoal
HCM222	SS IV Pole W/4 Hook Top, No Lose Thumb Knob, 5 Leg 16"dia. 22 lb. Base W/3" Casters, Charcoal
HCM223	SS IV Pole W/6 Hook Top, No Lose Thumb Knob, 5 Leg 16"dia. 22 lb. Base W/3" Casters, Charcoal

# IV Poles and Accessories



## HCM270 Series

### 5 Leg Spider IV Poles

This 5 Leg spider base IV Pole was designed to maximize stability. Base is 22" diameter and has a low center of gravity with enhanced weight for added stability. 3" scuff resistant casters are standard for ease of mobility. Height adjusts from 50"-92" with no-lose thumb knob. Available in chrome finish or durable stainless steel finish. Different top options are available to meet your needs. Designed to hold several infusion pumps. Standard epoxy coated base color is charcoal.

Colored bases are available at no extra charge to help each department maintain their assets. Simply add a – to the part number and follow with the first letter of the color.

Standard colors are:

<span style="background-color: black; width: 10px; height: 10px; display: inline-block;"></span>	CHARCOAL	<span style="background-color: red; width: 10px; height: 10px; display: inline-block;"></span>	R: RED	<span style="background-color: blue; width: 10px; height: 10px; display: inline-block;"></span>	B: BLUE
<span style="background-color: yellow; width: 10px; height: 10px; display: inline-block;"></span>	Y: YELLOW	<span style="background-color: darkgreen; width: 10px; height: 10px; display: inline-block;"></span>	G: GREEN		

#### Chrome Poles

HCM272	Chrome IV Pole W/Thumb Knob, 2 Hook Top, 5 Leg Spider Base W/3" Casters, Charcoal
--------	---

HCM273	Chrome IV Pole W/Thumb Knob, 4 Hook Top, 5 Leg Spider Base W/3" Casters, Charcoal
--------	---

#### Stainless Steel Poles

HCM270	SS IV Pole W/Thumb Knob, 2 Hook Top, 5 Leg Spider Base W/3" Casters, Charcoal
--------	---

HCM271-RT	SS IV Pole W/Thumb Knob, 8 Hook Rake-Top, 5 Leg Spider Base W/3" Casters, Charcoal
-----------	--

HCM271	SS IV Pole W/Thumb Knob, 4 Hook Top, 5 Leg Spider Base W/3" Casters, Charcoal
--------	---

HCM271-6	SS IV Pole W/Thumb Knob, 6 Hook Top, 5 Leg Spider Base W/3" Casters, Charcoal
----------	---



## HCM270 Series

### 6 Leg Spider IV Poles

This 6 Leg spider base IV Pole is our most sturdy IV Pole. Base is 22" diameter and has a low center of gravity with enhanced weight for added stability. 3" scuff resistant casters are standard for ease of mobility. Height adjusts from 50"-92" with no-lose thumb knob. Available in chrome finish or durable stainless steel finish. Different top options are available to meet your needs. Designed to hold multiple infusion pumps. Standard epoxy coated base color is charcoal.

Colored bases are available at no extra charge to help each department maintain their assets. Simply add a – to the part number and follow with the first letter of the color.

Standard colors are:

<span style="background-color: black; width: 10px; height: 10px; display: inline-block;"></span>	CHARCOAL	<span style="background-color: red; width: 10px; height: 10px; display: inline-block;"></span>	R: RED	<span style="background-color: blue; width: 10px; height: 10px; display: inline-block;"></span>	B: BLUE
<span style="background-color: yellow; width: 10px; height: 10px; display: inline-block;"></span>	Y: YELLOW	<span style="background-color: darkgreen; width: 10px; height: 10px; display: inline-block;"></span>	G: GREEN		

#### Chrome Poles

HCM277	Chrome IV Pole W/Thumb Knob, 2 Hook Top, 6 Leg Spider Base W/3" Casters, Charcoal
--------	---

HCM278	Chrome IV Pole W/Thumb Knob, 4 Hook Top, 6 Leg Spider Base W/3" Casters, Charcoal
--------	---

#### Stainless Steel Poles

HCM275	SS IV Pole W/Thumb Knob, 2 Hook Top, 6 Leg Spider Base W/3" Casters, Charcoal
--------	---

HCM276-RT	SS IV Pole W/Thumb Knob, 8 Hook Rake-Top, 6 Leg Spider Base W/3" Casters, Charcoal
-----------	--

HCM276	SS IV Pole W/Thumb Knob, 4 Hook Top, 6 Leg Spider Base W/3" Casters, Charcoal
--------	---

HCM279	SS IV Pole W/Thumb Knob, 6 Hook Top, 6 Leg Spider Base W/3" Casters, Charcoal
--------	---

# IV Poles and Accessories

## IV Pole Buying Guide

### Pole Selections

#### Standard



#### Chrome Pole:

For use in non-acute areas.

#### Premium



#### Stainless Steel Pole:

For use in OR and our most durable pole option

### Caster Selections

#### Standard



#### 2" Ball Bearing Casters:

Quality Caster for Lighter Weight Bases.

#### Premium



#### 3" Ball Bearing Casters:

For heavy duty bases. Better elevator transitions.

### Base Selections

#### Standard



#### 5 Leg Economy Base:

24" diameter

For lighter use areas.



#### 6 Leg Economy Base:

24" diameter

Better stability for use where only standard IV is administered.

#### Premium



#### 5 Leg Weighted Base:

16" diameter

Available in five standard colors. For infusion pumps when paired with SS Pole.



#### 5 Leg Spider Low Center of Gravity Base:

22" diameter

Available in five standard colors. For infusion pumps when paired with SS Pole.



#### 6 Leg Spider Low Center of Gravity Base:

22" diameter

Available in five standard colors. For infusion pumps when paired with SS Pole. Most stable pole and best elevator transition.

### Top Selections

#### Standard



#### 2-Hook Ram's Horn Top:

For light use areas.



#### 4-Hook Ram's Horn Top:

For standard applications.

#### Premium



#### 6-Hook Ram's Horn Top:

For use in specialty applications where several different IV Bags are needed.



#### 8-Hook Rake Top:

For use in specialty applications where several different IV Bags are needed.

# IV Poles and Accessories

## IV Pole Buying Guide

Item #	POLE SELECTION		CASTER SELECTION		BASE SELECTION				TOP SELECTION			
	Standard	Premium	Standard	Premium	Standard		Premium		Standard		Premium	
	Chrome Pole: For use in non-acute areas.	Stainless Steel Pole: For use in OR and our most durable pole option.	2" Ball Bearing Casters: Quality Caster for Lighter Weight Bases.	3" Ball Bearing Casters: For heavy duty bases. Better elevator transitions.	5 Leg Economy Base: 24" diameter For lighter use areas.	6 Leg Economy Base: 24" diameter Better stability for use where only standard IV is administered	5 Leg Weighted Base: 16" diameter Available in five standard colors. For infusion pumps when paired with SS Pole.	5 Leg Spider low center of Gravity Base: 22" diameter Available in five standard colors For infusion pumps when paired with SS Pole. Most stable pole and best elevator transition.	2-Hook Ram's Horn Top: For light use areas.	4-Hook Ram's Horn Top: For standard applications.	6-Hook Ram's Horn Top: For use in specialty applications where several different IV Bags are needed.	8-Hook Rake Top: For use in specialty applications where several different IV Bags are needed.
HCM201	X			X		X				X		
HCM202	X			X		X					X	
HCM210		X	X		X					X		
HCM210-RT		X	X		X							X
HCM211		X	X		X						X	
HCM212		X	X		X							X
HCM220	X			X			X			X		
HCM221		X		X			X			X		
HCM222		X		X			X				X	
HCM222-RT		X		X			X					X
HCM223		X		X			X					X
HCM234	X			X		X				X		
HCM235	X			X		X					X	
HCM237		X		X		X				X		
HCM237-RT		X		X		X						X
HCM238		X		X		X					X	
HCM239		X		X		X						X
HCM270	X			X				X		X		
HCM271		X		X				X			X	
HCM271-RT		X		X				X				X
HCM271-6		X		X				X				X
HCM272	X			X				X		X		
HCM273	X			X				X			X	
HCM275		X		X					X	X		
HCM276		X		X					X		X	
HCM276-RT		X		X					X			X
HCM277	X			X					X	X		
HCM278	X			X					X		X	
HCM279		X		X					X			X

# IV Poles and Accessories



## HCM290 Series

### 6 Leg Stainless Steel Spider Base IV Poles

This IV Pole is completely made of stainless steel to inhibit any type of corrosion or rust. 24" Diameter Base is heavy duty weighted stainless steel with low center of gravity. 1" stainless steel pole adjust from 50"-92" with no-lose thumb knob. Available with any of our 4 top options. Premium sealed ball bearing casters make this pole easy to maneuver and very stable.

- › Premium casters
- › Fully welded stainless steel construction



## HCM280 Series

### Double IV Pole with Steering Wheel and Tray

This IV Pole is designed for areas where several infusion pumps are needed such as ICU. Double poles make more spots to hang pumps and the two four hook tops make it easier to keep IV lines organized. Height adjusts from 85"-110". Our 6 Leg spider base with 3" ball bearing casters make this a very stable and easy to move pole. Full circle steering wheel makes it easy to maneuver from any side. Complete pole construction is 304 grade stainless steel. Removable plastic tray makes a great place for storage and a spot to set things during bag changes and setup.



## HCM290 Series

### 5 Leg Foot Control IV Pole

Fully welded stainless steel constructed IV Pole with 22" diameter 5 Leg spider base. Foot control height adjustment leaves hands free to adjust pole, IV Tubing and bags. Our 2" swivel **sayno2rust** casters make moving easy. Bar to protect foot pedal from getting bumped. 4 Hook ram's horn style top. 1.25" diameter pole adjusts in height from 52"-94". Easy to raise, simply by lifting pole to desired height. Ships fully assembled-no further assembly required.



## HCM200 Series

### Rigid IV Pole

An IV Pole that attaches directly to the OR Table rail. Designed to hold IV near the point of use for convenience as well as to save space in the OR. Constructed of durable stainless steel. Hook design allows bags to be changed rapidly. Pole measures 27 1/2" high. 1/2" diameter mounting post attaches to OR table rail with standard clark socket (*sold separately*).

HCM295-2	SS IV Pole W/Thumb Knob, 2 Hook Top, 6 Leg SS Spider Base W/3" Ball Bearing Casters
HCM295	SS IV Pole W/Thumb Knob, 4 Hook Top, 6 Leg SS Spider Base W/3" Ball Bearing Casters
HCM295-6	SS IV Pole W/Thumb Knob, 6 Hook Top, 6 Leg SS Spider Base W/3" Ball Bearing Casters
HCM295-RT	SS IV Pole W/Thumb Knob, 8 Hook Rake-Top, 6 Leg SS Spider Base W/3" Ball Bearing Casters

HCM286	Double IV Pole with Steering Wheel and Tray
--------	---

HCM290	5 Leg Foot Control IV Pole
--------	----------------------------

Upgraded casters included.  
**sayno2rust**®

HCM215	Rigid IV Pole
--------	---------------

Attaches to Clark Socket (*sold separately*)

**SEE PAGE 34**

# IV Poles and Accessories



## HCM200 Series Lift Assist IV Poles

An IV Pole designed for the nurse! Eliminates lift bags up high and straining backs and shoulders. Simply load the pole up with bags and turn the crank to raise. Adjusts in height from 56"-96" just by turning the crank. Large 6 leg base adds the needed stability. 3" scuff resistant casters provide unmatched maneuverability. Base post is 2" square stainless steel tube with upper extension tube of 1 1/4" stainless steel tube. No electric motors or hydraulic pumps to replace. Available with 4-hook, 8-hook or 28-hook top.

HCM251	4 Hook Stainless Steel Lift Assist IV Pole
HCM252	8 Hook Stainless Steel Lift Assist IV Pole
HCM253	28 Hook Stainless Steel Lift Assist IV Pole



## HCM200 Series Uni-Pole

This IV pole attaches to any round or square tube up to 1 1/4". Universal clamp design allows 180 degree rotation to set the pole wherever you need it. Adjust in height with integral clamp. Constructed of stainless steel for ease of cleaning. No additional clamps needed.

HCM216	Uni-Pole (adjust from 18"-36")
HCM216-4	Uni-Pole with 4 hook top (adjust from 18"-36")



## HCM200 Series IV Towbar

Constructed of stainless steel for transport of IV Poles, this towbar can transport two IV Poles at the same time. It can also drop into the hole on most stretchers and wheelchairs. Adjusts from 20"-30".

HCM240	Telescoping IV Towbar
HCM245	Universal Clamp with Hole for Attaching Towbar



## HCM200 Series Universal Clamp for IV Pole Accessories

This aluminum clamp attaches most of our IV pole accessories to the IV Pole. It will also clamp to tables, wheelchairs, carts or any other round or square tube up to 1 1/2" in diameter.

HCM231	Universal Clamp for IV Pole Accessories
--------	--

# IV Poles and Accessories



## HCM200 Series

### Universal Clamp for Lift Assist IV Pole

HCM257

Universal Clamp for  
Lift-Assist IV Pole

This clamp attaches most of our IV pole accessories to the lift assist IV Pole. Works much the same as our HCM231 but for the Lift Assist IV Pole. Makes any accessory work with the Lift Assist Pole.



## HCM200 Series

### Power Strip for IV Pole

HCM232

Power Strip for IV  
Pole

Fit out with 6 Hospital Grade Receptacles and an integral clamp, this is a great IV Pole Power Strip. Measures 5.65" H x 5.28" W x 3.27" D. Meets UL-1363-A for use in patient areas. Fits IV pole diameter 7/8"-1 1/2". 15 foot hospital grade power cord attached.



## HCM200 Series

### Add-A-Tray for IV Pole

HCM260

Add-A-Tray  
for IV Pole

This Stainless Steel Tray measures 13" x 9" and attaches to the IV Pole with our Universal Clamp HCM231(*sold separately*). Excellent spot to put small instruments or items during set up. Stainless steel construction makes it easy to clean. With our universal clamp can easily mount to tables, mayo stands or carts.

Mounts with Universal Clamp  
(*sold separately*)

[SEE PAGE 53](#)



## HCM200 Series

### Infusion Pump Frame for IV Pole

HCM259

Infusion Pump Frame  
for IV Pole

This frame attaches to an IV Pole with two universal clamps (*sold separately*). Designed to hold up to 4 infusion pumps. 1" tubular frame in chrome construction. Measures 17 3/4" x 27 3/4".

Mounts with 2 Universal Clamps  
(*sold separately*)

[SEE PAGE 53](#)



## HCM200 Series

### Oxygen Tank Holder for IV Pole

HCM229

Oxygen Tank Holder  
for IV Poles

Heavy duty stainless steel construction for holding Oxygen tanks on IV poles or anywhere the universal clamp can attach to. Attaches with 1 universal clamp (*sold separately*). Accommodates most styles of oxygen tanks up to 4 7/8" diameter.

Mounts with Universal Clamp  
(*sold separately*)

[SEE PAGE 53](#)

# IV Poles and Accessories



## HCM200 Series

### Wire Basket for IV Pole

HCM217

Wire Basket  
for IV pole

Mounts with Universal Clamp  
(sold separately)

[SEE PAGE 53](#)



## HCM200 Series

### Stat Hanger for IV Pole

HCM200

Stat Hanger for  
IV Pole

Mounts with Universal Clamp  
(sold separately)

[SEE PAGE 53](#)



## HCM200 Series

### Urinal Bag Hook

HCM258

Urinal Bag Hook

Stainless steel hook designed to hold the Foley catheter bag lower on the IV Pole. Clamps directly to 1" pole. No other clamps needed.



## HCM200 Series

### Steering Wheel and Patient Tray for IV Pole

HCM230

Steering Wheel

HCM243

Patient Tray

16" diameter epoxy coated steering wheel attaches to 1" IV Pole. Makes it easy to steer pole from any side. Gives patient added stability in transport. Large diameter makes it easy for patient and caregiver to hold onto pole. Can be adjusted to any height. Easy to clean, plastic patient tray snaps easily into triangle on steering wheel.



## HCM200 Series

### Hand Sanitizer Holder for IV Pole

HCM261

Hand Sanitizer Holder  
for IV Pole

Keep the hand sanitizer handy on your IV Pole or wheelchair. Attaches to pretty much any surface without any other clamps needed. Easily adjusts to various sizes of hand sanitizer bottles. Durable plastic construction.

# Index

## A

Accessories for OR Tables	29-38
Accessory and Stirrup Carry Cart	38
Add-A-Rail	45
Additional OR Equipment	19-22
Anesthesia Cart	9
Anesthesia Screen	
Adjustable	35
Malleable	35
Rigid	35
Shield	35
Arm and Hand Surgery Tables	39-46
Arm and Hand Tables	
Carbon Fiber	
Rectangle	40
Under Pad	43
Carpal	40
Free Access	40
Hourglass Carbon Fiber	41
Hourglass Phenolic	41
Hourglass Under Pad	42
Minor Procedure	40
Mobile Base	
Rectangle	44
Hourglass	44
Rectangle Carbon Fiber	41
Rectangle Phenolic	41
Rectangle Under Pad	42
Stainless Steel Under Pad	43
Armboards	
Droplatch	
Bariatric	31
Standard	30
Pads	33
Radiolucent Height Adjustable	30
Rail Locking Droplatch	32
Under Pad	31
Extra Wide Under Pad	31
Wall Mount	33

## C

Candy Cane Stirrup	36
Carry Cart	
Step Stool	21
Stirrup	21
Clark Socket	34
Elite	34

## D

Double IV Pole	52
Double Picket Fence Leg Holder	37
Drop Leaf Instrument Table	11

## E

Elbow Arthroscopy	33
Elite Clark Socket	34
Equalizing Device	45
EZ-Tranz Gold Patient Transfer Boards	26
Replacement Covers	26
Wall Rack	27

## G

Glove Boxes	28
Gowning Benches	22

## H

Hampers	
Round	25
Square	25
Triangle	25
Horizontal Traction Towers	45
Hose Tree	35

## I

Insta-Rail	38
Instrument Tables	
Drop Leaf	11
Height Adjustable	11
Nested	7
Two Tier Back Table	11
with Drawers	10
with H-brace	7
with Shelf	7
IV Poles	

5 Leg Space-Saving	48
5 Leg Foot Control	52

Economy	
5 Leg	48
6 Leg	48

Lift Assist	53
-------------	----

Spider	
--------	--

5 Leg	49
6 Leg	49
6 Leg Stainless Steel	52

## IV Pole Accessories

Add-A-Tray	54
Power Strip	54
Rigid	52
Hand Sanitizer Holder	55
Infusion Pump Frame	54
Oxygen Tank Holder	54
Stat Hanger	55
SteeringWheel and Patient Tray	55
Urinal Bag Hook	55
Wire Basket	55
Universal Clamp	53
IV Poles and Accessories	47-55
IV Pole Buying Guide	50-51
IV Towbar	53

## K

Kick Bucket	20
-------------	----

## L

Lift Assist Tables	16
--------------------	----

## M

Mayo Stands	13-18
Foot Control	14
Knob Control	15
Thumb Control	15
Aluminum Base Mayo Stand	17
Stainless Steel Base	
Mayo Stand	17
Mayo Tray Attachment	17, 38

## P

Patient Transfer Roller Boards	26
Replacement Covers	26
Wall Rack	27
Picket Fence Leg Holder	37

## R

Replacement Pads	45
Rigid IV Pole	52

## S

Side Rail Extension	38
Side Rail Socket	
Aluminum	34
Stainless Steel	34

## Single Loop Leg Prep with Pad

with Pad	36
----------	----

## Slider Style Patient Transfer Boards

Boards	26
--------	----

## Wall Rack

Wall Rack	27
-----------	----

## Solution Stands

Double Bowl	20
-------------	----

## Single Bowl

Single Bowl	20
-------------	----

## Spider Base IV Poles

6 Leg	49
-------	----

## Sponge Receptacle

Sponge Receptacle	20
-------------------	----

## Stainless Steel Tables

Stainless Steel Tables	5-12
------------------------	------

## Stand Alone Leg Prep

Stand Alone Leg Prep	36
----------------------	----

## Step Stools

Chrome	21
--------	----

Stacking Interlocking	21
-----------------------	----

## Sterile Wrapping Table

Sterile Wrapping Table	10
------------------------	----

## Stools

Physicians	22
------------	----

Surgical	22
----------	----

## Storage, Transfer and Transport

Storage, Transfer and Transport	23-28
---------------------------------	-------

## Surgical Stools

Surgical Stools	22
-----------------	----

## T

## Tables for over OR Tables

Tables for over OR Tables	11
---------------------------	----

## TKR Support

TKR Support	37
-------------	----

## Two-tier Back Table

Two-tier Back Table	11
---------------------	----

## U

## Uni-Pole

Uni-Pole	53
----------	----

## Universal Clamp

for Lift Assist IV Pole	54
-------------------------	----

## Universal Rail Clamp

Universal Rail Clamp	34
----------------------	----

## Utility Cart with Baskets

Utility Cart with Baskets	24
---------------------------	----

## Utility Carts

Utility Carts	24
---------------	----

## Utility Tables

Utility Tables	8
----------------	---

with 4 sided Guardrail	9
------------------------	---

Large Utility/Prep	9-10
--------------------	------

## W

## Wall Shelves

Wall Shelves	28
--------------	----

## Width Extender

Width Extender	32
----------------	----

Bariatric Width Extender	32
--------------------------	----

## Work Tables

Work Tables	12
-------------	----